MS Schippers

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
Elektrisch schema

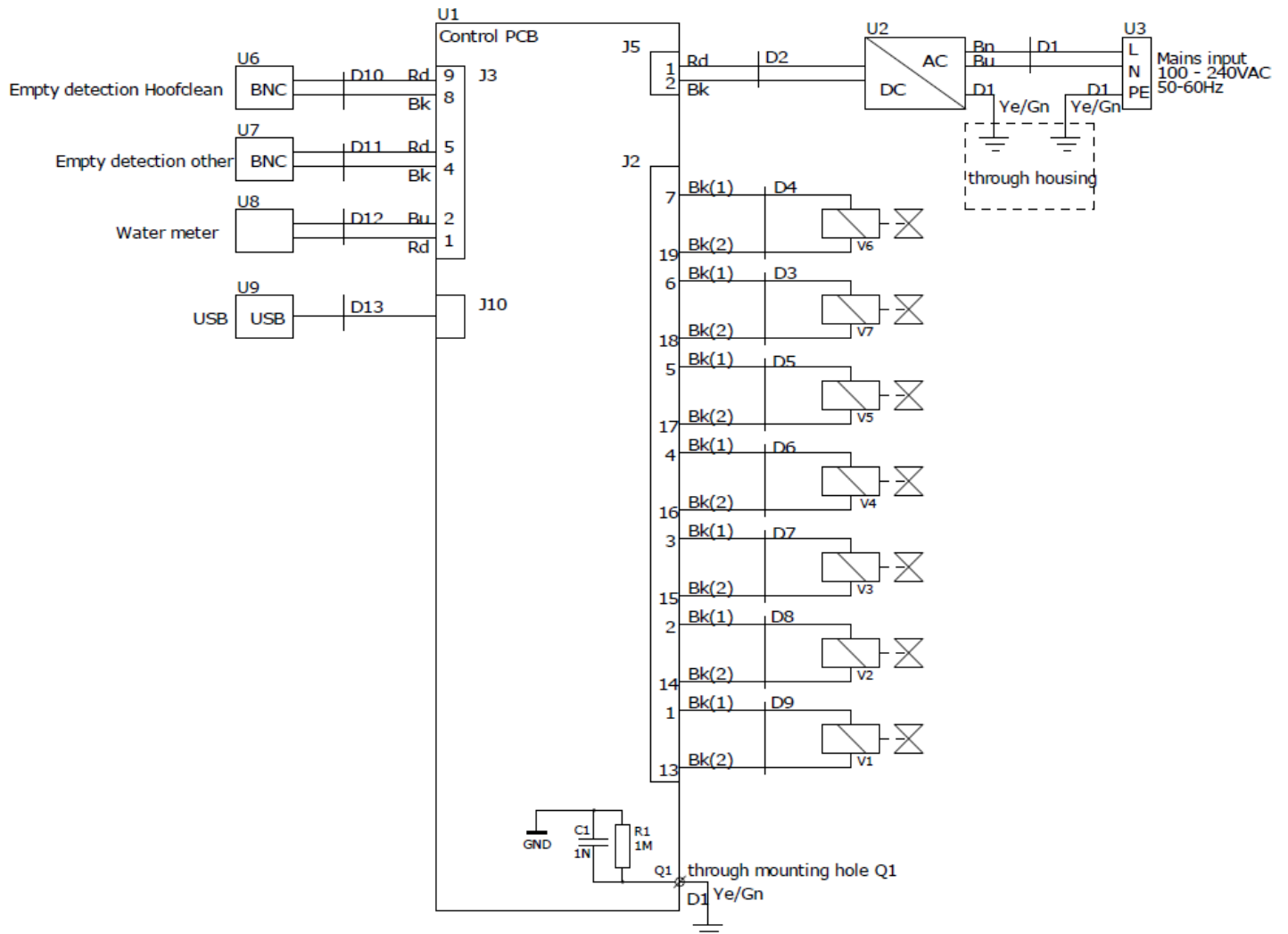
Contents

1. Quick start setup, milking parlor.....	57
2. Quick start setting, milking robot.....	58
3. General information.....	59
3.1 Liability & Conditions.....	59
3.2 General safety instructions.....	59
4. Introducing the AutoHoofClean pump.....	61
4.1 AutoHoofClean pump parts.....	62
5. Correct setup of the AutoHoofClean pump, hoses and product packaging.....	63
6. Commissioning and use.....	65
6.1 Installing the AutoHoofClean pump.....	65
6.2 Mounting the filter housing.....	67
6.3 Connecting the water supply to the AutoHoofClean pump.....	68
6.4 Installing the suction hoses.....	69
6.5 Installation of the vacancy detector and hose pillar.....	72
6.6 Connecting the mats to the AutoHoofClean pump.....	75
6.8 Additional points of interest.....	76
6.9 Connecting power supply.....	77
7. Display and system settings.....	78
7.1 Setting language and time.....	80
8. Startup and testing.....	81
8.1 Select stable type.....	83
8.2 Milking robot settings.....	84
8.3 Troubleshooting – milking robot.....	85
8.4 Setting Passages.....	86
8.5 Set other product.....	87
9. Milking parlor settings.....	88
9.1 Starting milking – setting the function.....	89
9.2 Start automatic milking.....	90
9.3 Starting manual milking.....	90
9.4 Setting up the Claw Problem.....	91
9.5 Setting Passages.....	92
9.6 Set other product.....	93
10. Meanings of icons on screen.....	94
11. Connect the AutoHoofClean pump via the HyCare app.....	95
12. MS HoofClean Care & other product empty notification.....	96
13. Read data via USB port.....	97


14. Screen installer.....	98
14.1 Setting Language/Time.....	98
14.2 Counter screen.....	99
14.3 Alarm level.....	100
15. Settings in the control panel (only applicable to milking parlors).....	101
16. View technical data via the “Info” menu.....	102
17. Maintenance of the Autohoofclean pump.....	103
18. Replacing the water filter (hardware).....	104
18.1 Replacing the water filter (program).....	105
Completion of the installation.....	105

1. Quick start setup, milking parlor


1. Select gear

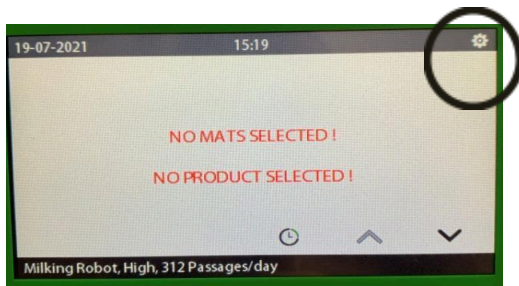
2. Choose "Claw Treatment"

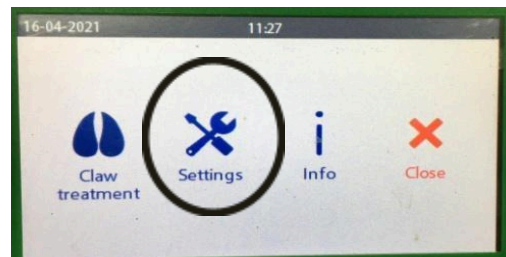
3. Unlock computer

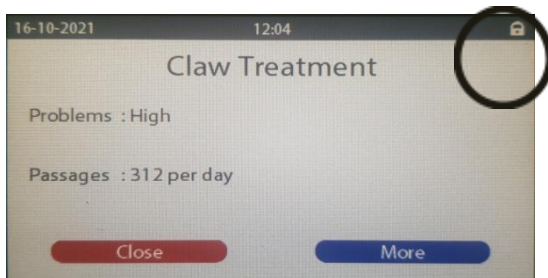4. Enter code

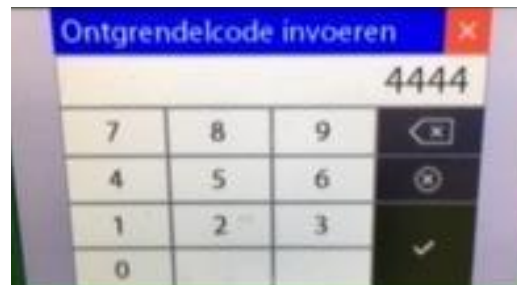
5. Select "barn Type"

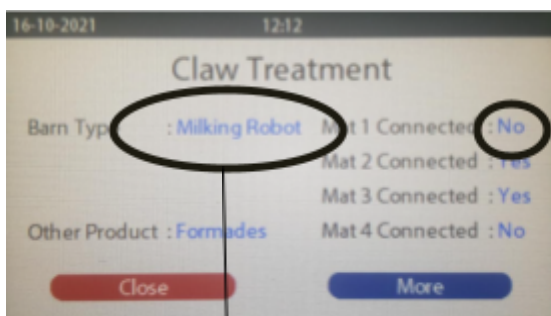
6. Select the mats and choose "more"


7. Select product

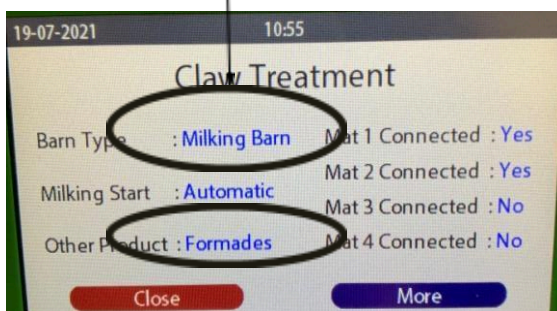
8. Enter the milking duration and start time

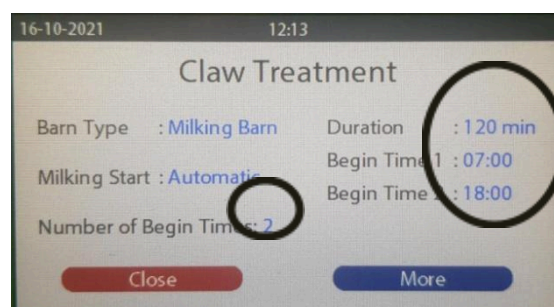9. Choose product

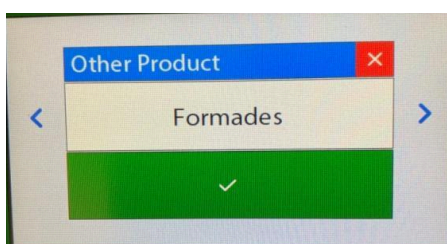
10. Select "Claw Treatment"

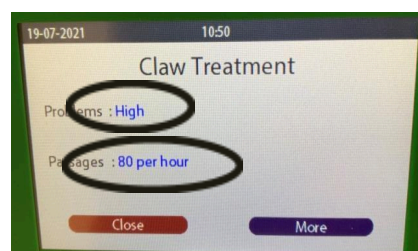
2. Quick start setting, milking robot


1. Start the AutoHoofClean pump

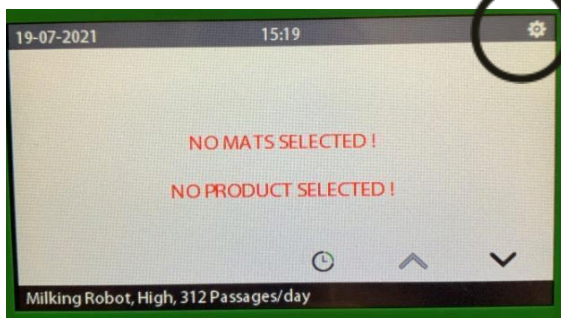2. Select "Claw Treatment"

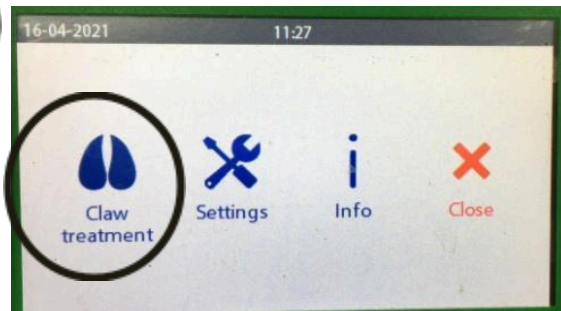
3. Unlock computer

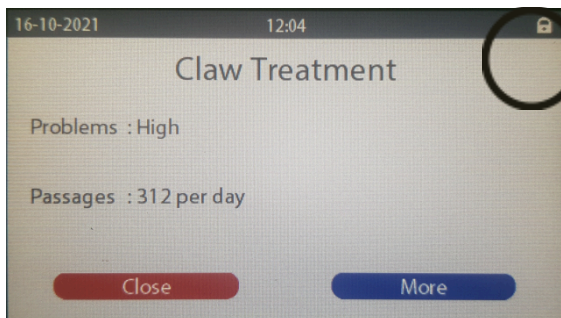4. Enter code

5. Choose product

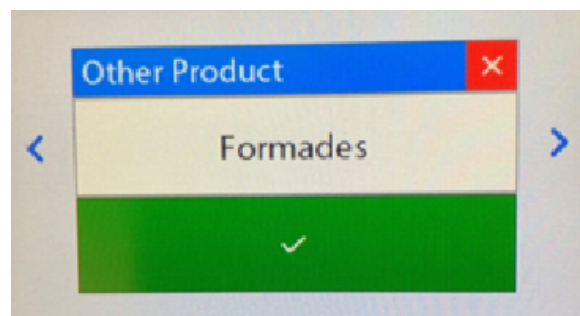
6. Set connected outputs to "Yes"

7. Press "more"

8. Set "Problems" to "High"

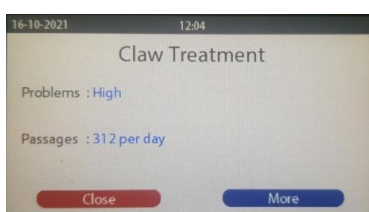3. General information

This manual is intended to support the use and maintenance of the AutoHoofClean pump. Do not use the AutoHoofClean pump until you have read and understood the entire manual. Contact The Schipper Group specialists if you have any questions about this user manual, its use, and/or maintenance. The following notes and explanations are essential for using the device.

3.1 Liability & Conditions

Only qualified personnel should operate the AutoHoofClean pump in accordance with all safety requirements, including all applicable laws, regulations, and standards.

This documentation has been carefully compiled, but the device described is undergoing continuous development. Therefore, it is possible that the documentation is not always complete or does not fully comply with the described performance data, standards, or other features. In the event of technical or typographical errors, we reserve the right to make changes to products already supplied at any time and without prior notice. MS Schippers is not responsible for the content of this document.

3.2 General safety instructions

The owner of the AutoHoofClean pump must ensure that the operating personnel are aware of the contents of this manual and that it is always available to them.

Keep the AutoHoofClean pump out of reach of children, pets, and unauthorized persons. The AutoHoofClean pump is designed for use in a cattle barn.

The user is responsible for operating the AutoHoofClean pump safely. The user must comply with all safety instructions before operating the AutoHoofClean pump. Wear appropriate personal protective equipment (clothing, goggles, dust mask, earmuffs, gloves).

NEVER allow untrained personnel to operate the AutoHoofClean pump.

Before use, hang the AutoHoofClean pump according to the instructions in this installation manual. It is not permitted to use the pump in any other position, for example lying down.

Always use grounded outlets with the following specifications: Voltage 100-240VAC. Frequency: 50/60 Hz. Power (max) 50W.

The AutoHoofClean pump must be connected to a water supply with a pressure between 2 and 6 bar. For higher pressures, a pressure regulator must be installed. At lower pressures, proper operation of the device cannot be guaranteed.

Use only clean, drinking-quality tap water. This appliance must be equipped with water that meets the quality requirements of the European Drinking Water Directive (EU 2020/2184). Furthermore, it is not permitted to use rainwater, surface water or water with additives (e.g. cleaning agents, antifreeze).

The appliance must be connected to the water supply in accordance with applicable national and European regulations (including NEN-EN 1717 and the European Drinking Water Directive 2020/2184). To protect the drinking water supply, a suitable backflow preventer must be installed where necessary to prevent backflow (INTERRUPTION). Consult your water supplier or installer for this.

4. Introducing the AutoHoofClean pump

Congratulations on the purchase of your AutoHoofClean pump.

The AHC pump is an innovative system that efficiently cleans and disinfects cows' hooves. It functions as a dosing system, automatically pumping MS HoofClean Care and a complementary product (such as MS HoofClean Protect or MS Formades) along with water into the AutoHoofClean mats.

Up to four mats can be connected simultaneously, with fully automatic filling.


Structure and operation

The AutoHoofClean pump is housed in a robust stainless steel housing and is operated via a user-friendly touchscreen. Thanks to integrated valves, the correct amounts of water, disinfectant, and cleaning agent are always delivered precisely to the mats.

This document provides you with a detailed description of how to use the device correctly.

4.1 AutoHoofClean pump parts

No.	Description:
1.	Handles
2.	Touchscreen
3.	Water filter
4.	Plug/chassis part
5.	USB port
6.	Geka coupling
7.	BnC Connector

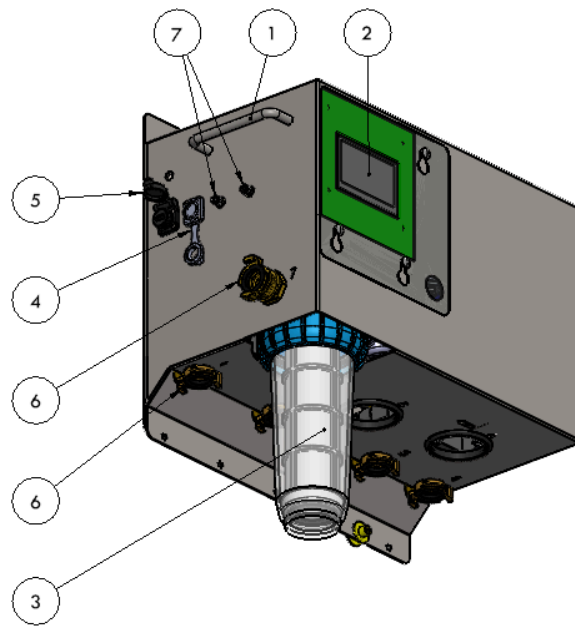4.2 Technical data AutoHoofClean pump

Device type	AutoHoofClean pump
System construction	Stainless steel housing
Complies with guidelines	CE
Dimensions	L=465mm W=269mm H=416mm
Weight	23 kg
Operating water pressure (min/max)	2-6 bar
Storage and operating temperature	5 - 35 °C

5. Correct setup of the AutoHoofClean pump, hoses and product packaging


For optimal operation of the AutoHoofClean pump, it is essential that the pump unit, product packaging and hoses are set up correctly.

Correct installation prevents malfunctions, such as siphoning, and ensures that the pump functions reliably.

⚠ Please note:

If the user chooses a different setup than described below, this is entirely at the user's own risk and responsibility.

Below you will find a schematic representation of the correct arrangement of the pump, packaging and hoses.

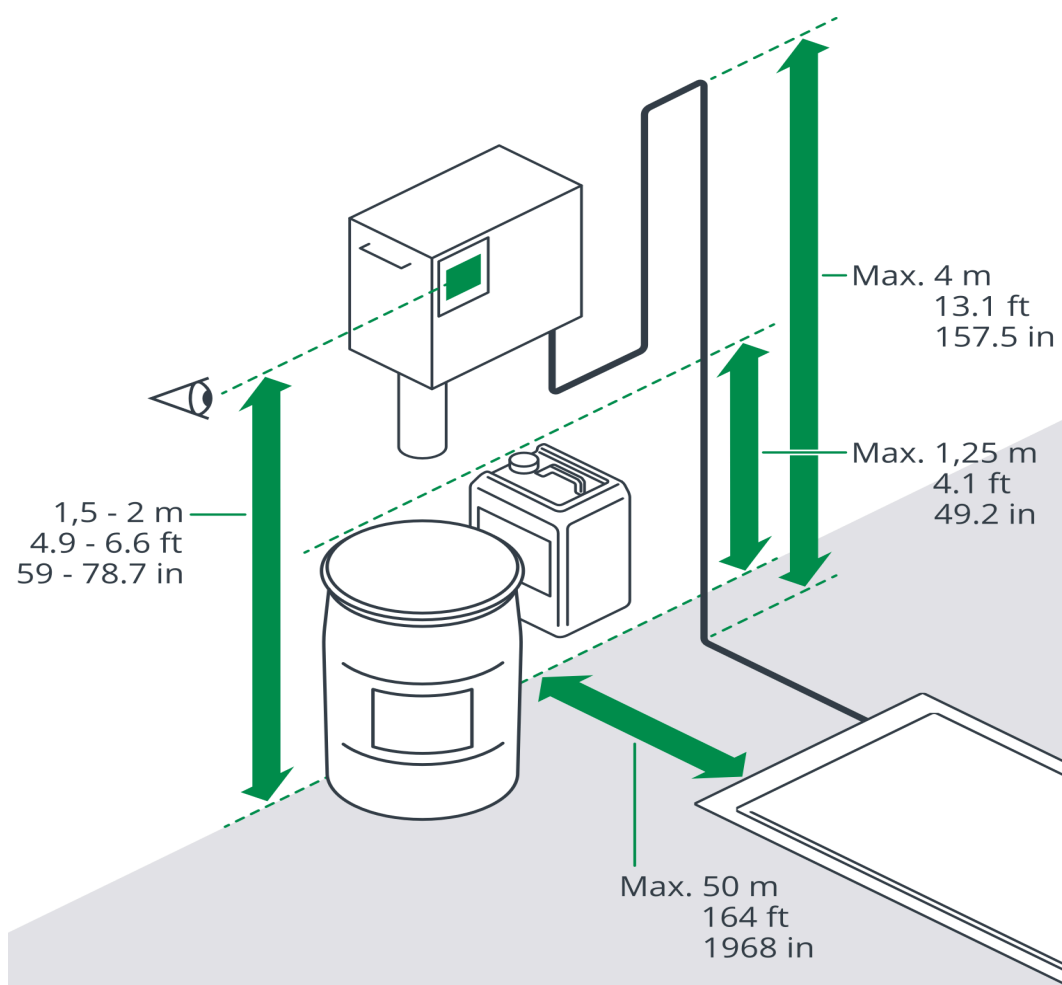1. Pump unit

- Placement at eye level
Mount the pump unit at eye level for easy operation and reading of the touchscreen.
- Maximum height difference with the mats
The height difference between the pump unit and the mat(s) may not exceed 2 metres.

✗ It is therefore not permitted to install the pump unit and/or product packaging in an attic, for example.

2. Product packaging (MS HoofClean Care & Protect/Formades)

- Level ratio with the mats
The cans should be approximately level with the mats.
 - Assembly of cans
Cans may be hung slightly off the ground (e.g. in a wall bracket), provided that the top of the can is no more than 1.25 metres above the mat level.
-

3. Snakes

- Maximum height difference
The height difference between the hose and the mat level may not exceed 4 metres.
 - Horizontal distance
The pump unit may be placed a maximum of 50 metres horizontally from the mat(s).
 - Total hose length
The total length of the hoses may not exceed 60 meters.
-

6. Commissioning and use

This chapter explains step-by-step how to prepare the AutoHoofClean pump for use and how the system works during normal use.

6.1 Installing the AutoHoofClean pump

Two bolts and plugs are supplied for correct and safe mounting of the AutoHoofClean pump.


⚠ Important:

Always mount the AutoHoofClean pump on a flat, vertical wall that is strong enough to support the weight of the pump.

When determining the mounting height, observe the conditions described in the chapter "Correct installation".

Step-by-step plan for installation:

1. Determine the mounting height
 - Choose a height where the pump is easily accessible and meets the requirements.
1. Mark the drill holes
 - Mark two drill holes horizontally at a centre distance of 350 mm.
2. Drill the holes
 - Use an 8mm drill bit to drill both holes.
3. Place the plugs
 - Hammer or push the supplied plugs into the drilled holes.
4. Tighten the bolts into the plugs
 - Screw the bolts into the plugs, but leave them sticking out a bit so you can easily hang the pump over them.
5. Hang up the AutoHoofClean pump
 - Place the AutoHoofClean-pump over the bolts and then tighten them firmly so that the pump is secure.

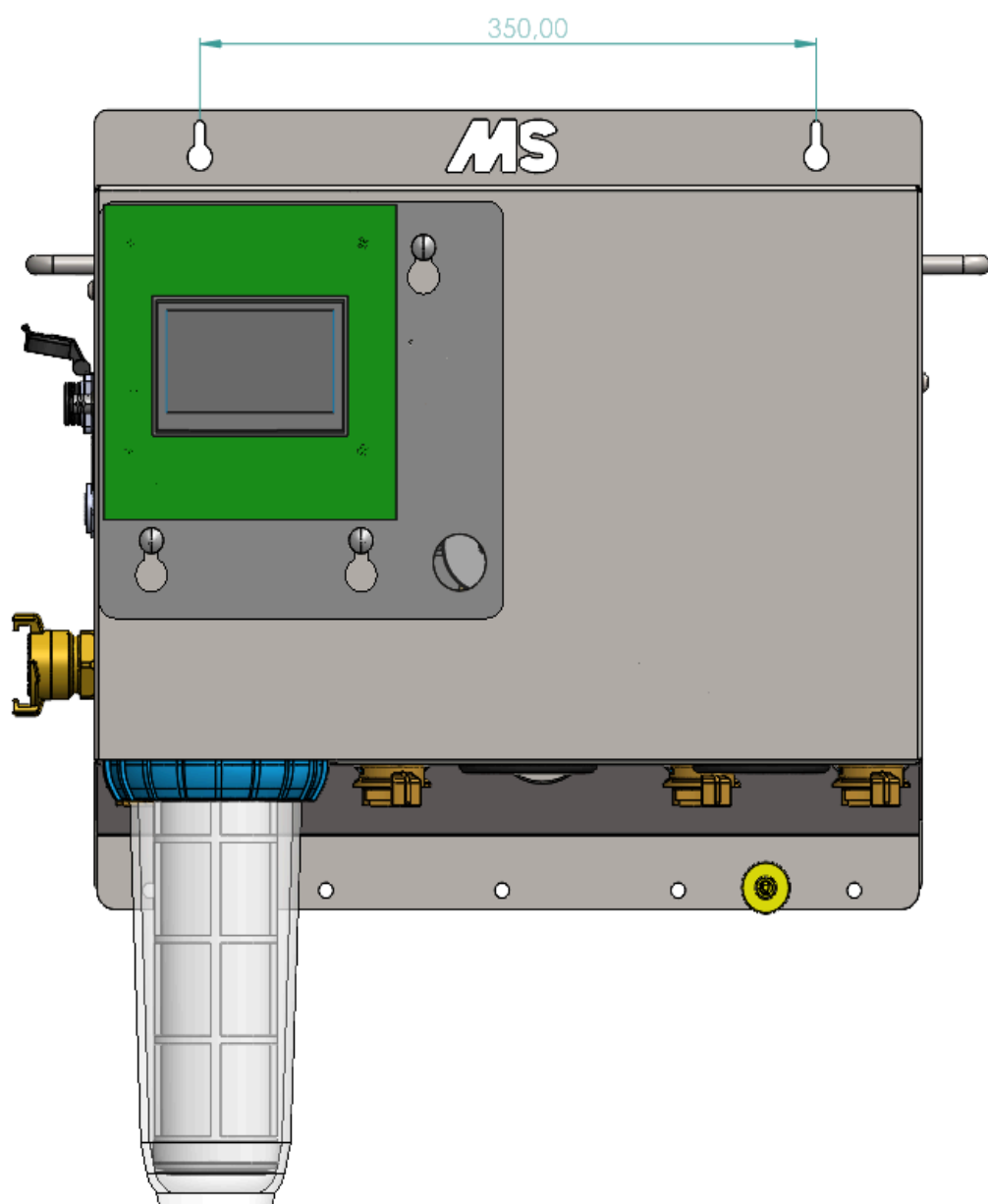
6.2 Mounting the filter housing

Before using the AutoHoofClean, the filter housing must be properly installed. Carefully follow the steps below:

Step-by-step plan:

1. Remove the packaging
 - Remove the filter cartridge from the plastic packaging, if you have not already done so.
2. Insert the filter cartridge
 - Place the cartridge upright in the filter housing.
3. Install the filter housing underneath the AutoHoofClean pump.
 - Screw the filter housing under the pump hand-tight.
4. Tighten the filter firmly
 - Use the supplied filter wrench to fully and firmly tighten the filter housing.

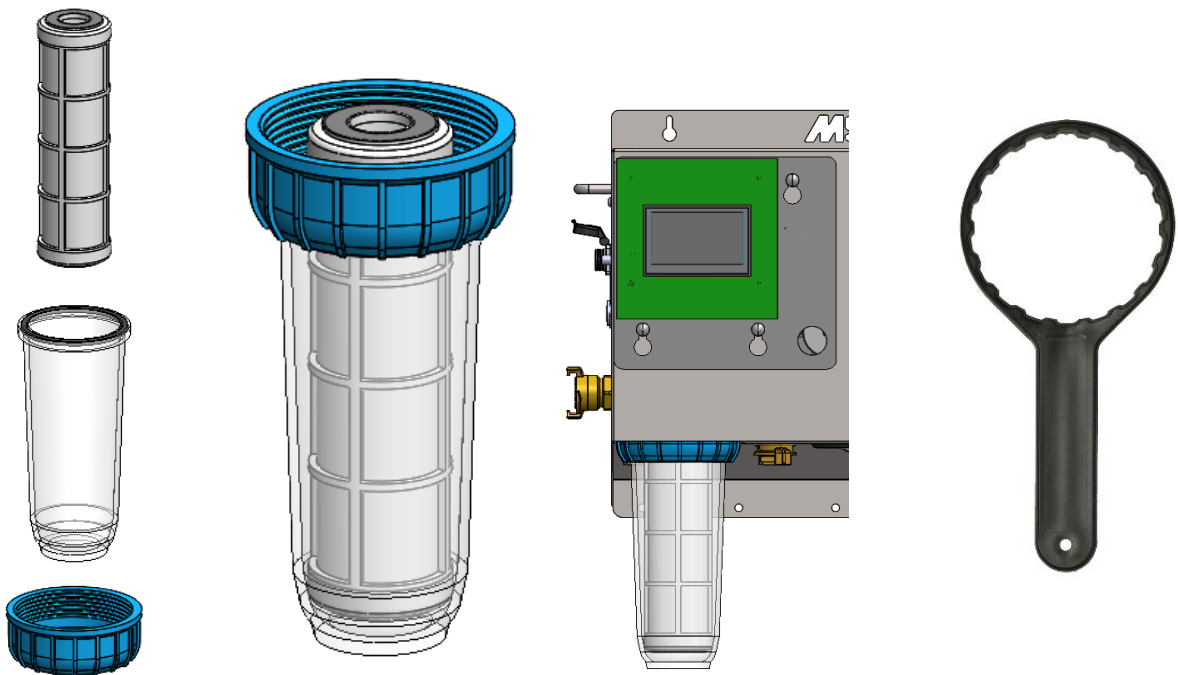
6.3 Connecting the water supply to the AutoHoofClean pump

For the AutoHoofClean pump to function properly, the water supply must be connected correctly. Follow the steps below carefully to ensure a leak-free and reliable connection.

Step-by-step plan:

1. Preparing the hose
 - Slide a hose clamp over the end of the water hose.
2. Attach to the GeKa coupling
 - Insert the water hose completely over the grommet of the GeKa coupling.
 - TIP: Soak the hose in warm water for a few minutes to soften it. This makes connecting easier.
 - Tighten the hose clamp securely to prevent leakage.
3. Connect to the AutoHoofClean pump
 - Connect the GeKa coupling to the water inlet of the AutoHoofClean pump.
 - Tighten the locking ring of the GeKa coupling firmly so that the coupling can no longer rotate.
3. Turn on the water supply
 - Turn on the water supply to allow water to flow to the pump.
4. Control
 - Check that there is at least 2 bar pressure in the system.
 - Inspect all connections for possible leaks.

✓ The water supply is now properly connected. You can proceed with setting up the pump or filling the mats.

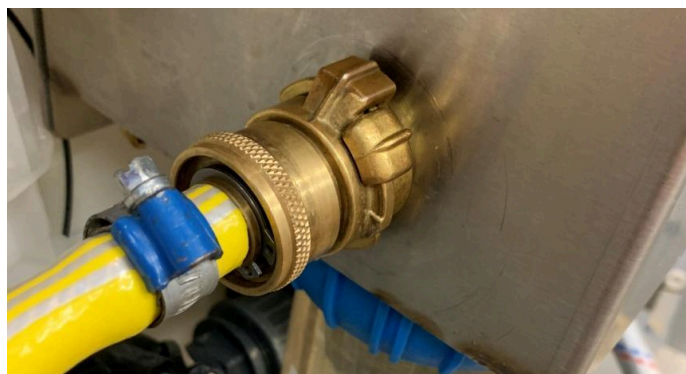6.4 Installing the suction hoses

To correctly dose the disinfectant/protective agent and the cleaning agent into the system, two suction hoses must be connected to the AutoHoofClean pump.

There are two connection points (holes) provided for this purpose at the bottom of the pump housing:


- Left hole (position 5 in the image below): for connection to MS HoofClean Care
- Right hole (position 6 in the image below): for connection to the second product

Important guidelines:

⚠ Use only the original suction hoses of the AutoHoofClean pump.

It is not permitted to extend the hoses or use alternative hoses.

- Standard length: 1.75 meters
- Available lengths: 2.5 meters and 3.5 meters
- Maximum permitted length: 3.5 meters

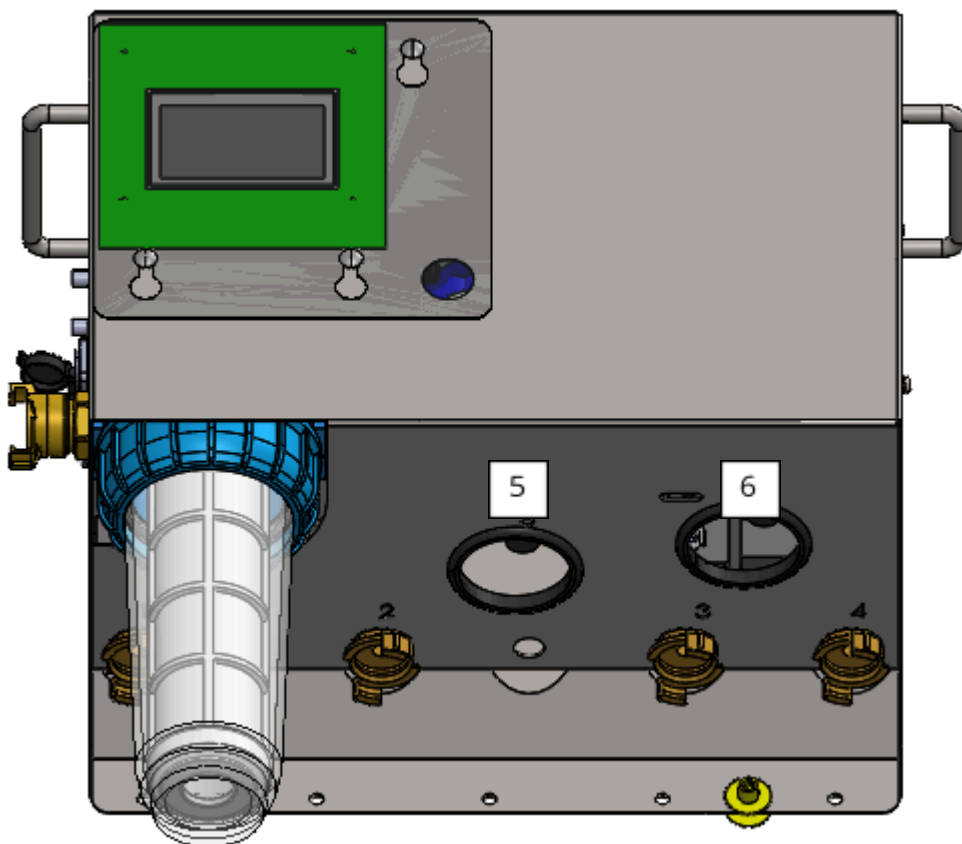When you look at the bottom of the AutoHoofClean pump, you will see the two connection points (pumps) to which the suction hoses need to be attached.


To properly use the AutoHoofClean pump, the supplied suction hoses must be connected to these pumps.

Step 1: Prepare the pump

- Unscrew the black nut from the pump.
- Remove the red protective cap from the pump.

This is a temporary barrier that was installed during transport.

- Keep the black nut safe, as you will need it in the next step to attach the hose.

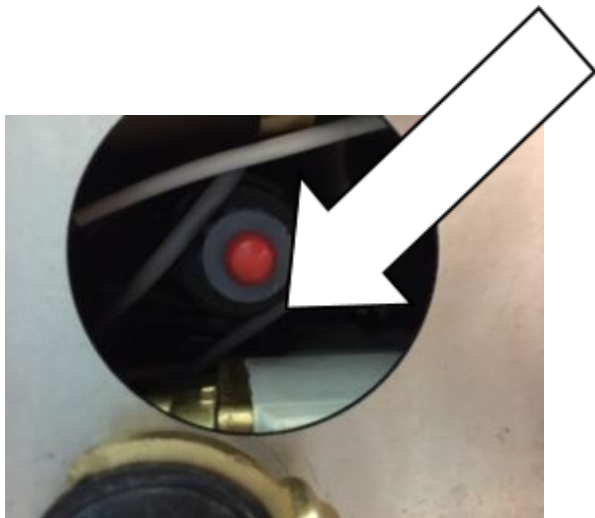Step 2: Preparing the suction hose

- Slide the black nut you removed in step 1 over the intake hose.

Step 3: Attaching the suction hose

- Insert the suction hose fully onto the pump until it stops moving.
Next, tighten the black nut firmly onto the pump to secure the hose.

Step 4: Check confirmation

- Check that both suction hoses are securely attached by gently pulling on each hose.
The hoses should remain firmly in place and not come loose.

6.5 Installation of the vacancy detector and hose pillar

To receive a timely warning when your barrels are (almost) empty, two empty status detectors and hose pillars are included.


This chapter explains how to correctly connect these to the transparent suction hoses and place them in the barrels.

Placement of the barrels

Step 1: Place the barrels in the previously described position under the pump.



Step 2: Drill a 20mm hole into the cap of your barrel.

- Then insert the cable from the vacancy detector together with the suction hose through this hole.

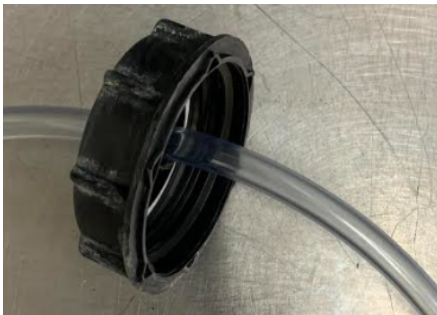
Step 3: Connecting the vacancy detector

- Connect the vacancy detector to the appropriate connectors on the AutoHoofClean pump:
 - Left connector (orange arrow): for MS HoofClean Care
 - Right connector (blue arrow): for the second product

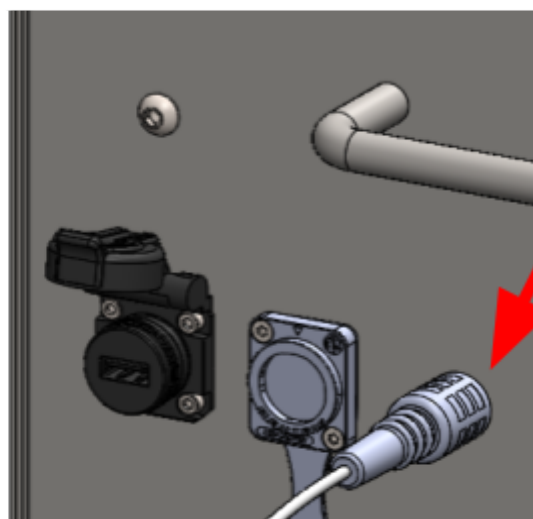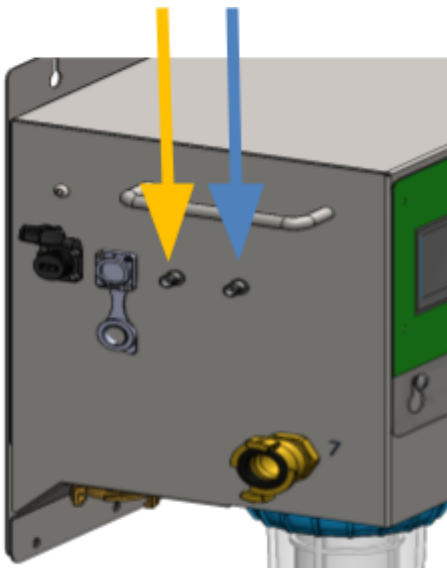Step 4: Mount the hose barb on the vacancy detector.

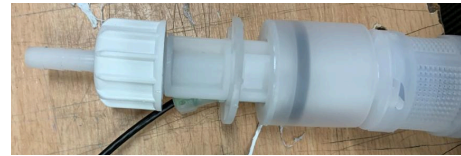Step 5: Slide the transparent suction hose over the hose pillar.

Step 6: Placing the vacancy detector

- Place the vacancy detector in the can and tighten the cap securely.
- Make sure the float is upright (vertically) with its bottom resting on the bottom of the barrel.


Step 7: Attaching the clamp

- Place the supplied clamp on the suction hose, as shown below. Make sure the clamp rests on the cap of the can.
- Please note that the distance between the bottom of the clamp and the top of the float is 24cm.

- Use the supplied M6x25 allen bolt, washer and M6 lock nut to secure the clamp.

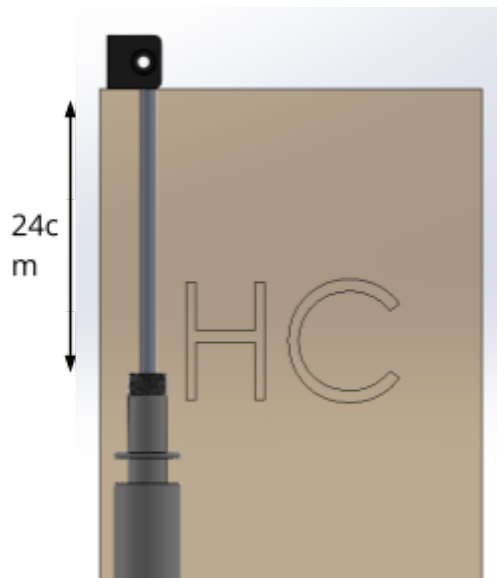
6.6 Connecting the mats to the AutoHoofClean pump


To supply your AutoHoofClean mats with product, they must be connected to the AutoHoofClean pump.

Up to four mats can be connected to the AutoHoofClean pump.

The GeKa couplings on the bottom of the pump are intended for the mats and are numbered from left to right: 1 to 4.

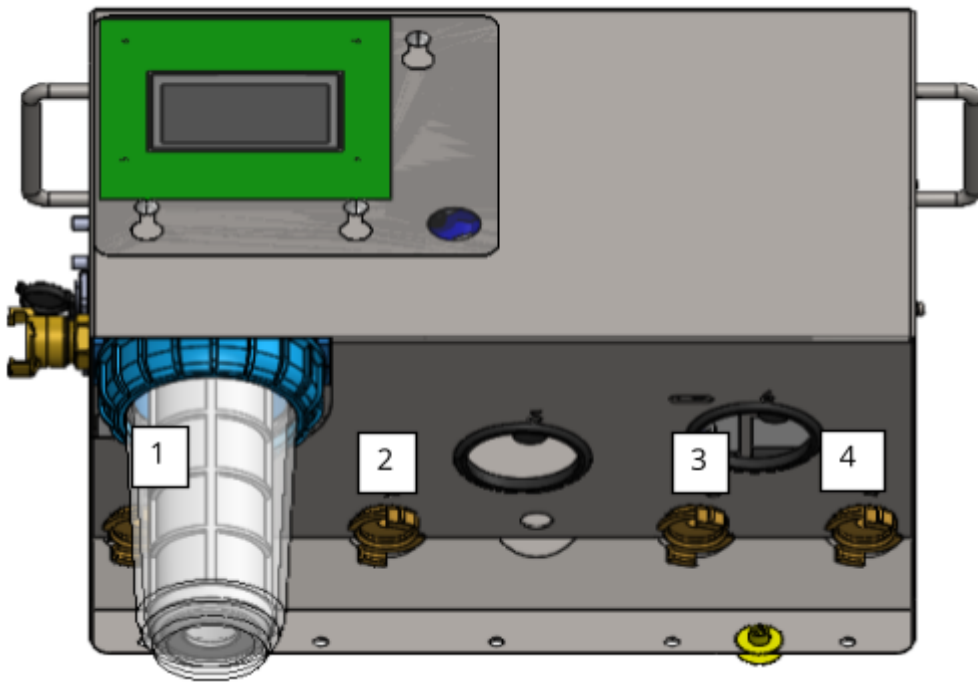Connect the GeKa couplings on the hoses leading to the mats to your chosen outlet. Then tighten the ring to ensure the coupling is leak-proof.

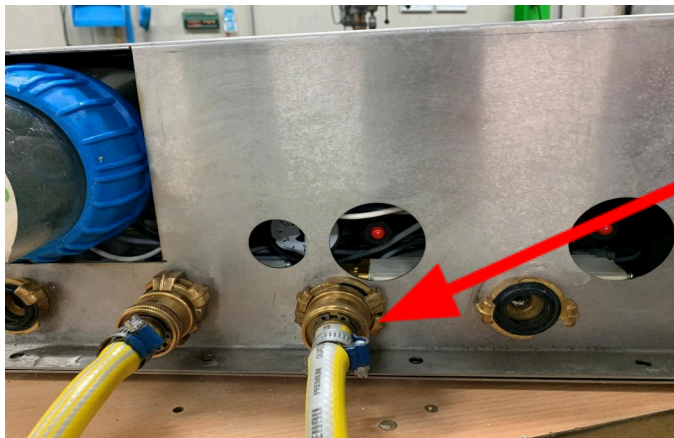
6.7 Placing the antenna on the MS AutoHoofClean

1. Connecting the antenna antenna

- Connect the antenna cable to the designated connection on the side of the MS AutoHoofClean pump (see image, marked in red) and tighten it by hand.
- Then mount the antenna (No. 1) on the bracket (No. 2) and connect the antenna cable (No. 5) to the antenna, tightening it by hand.

Please note: ensure that everything is mounted in the correct order; the order is shown below.

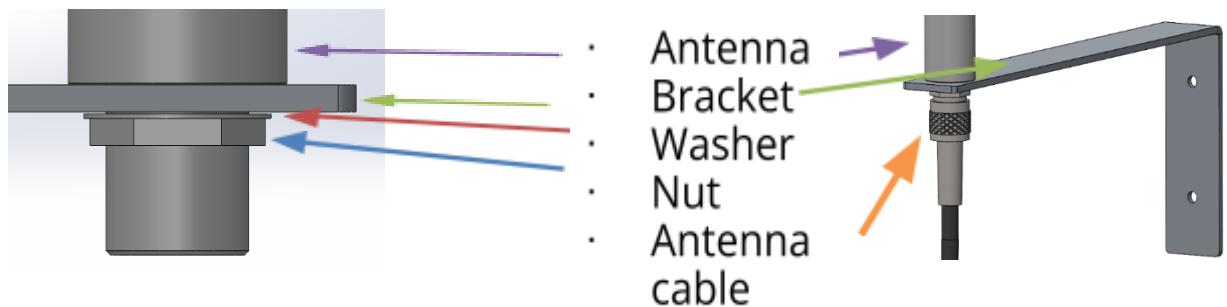2. One-time startup

- Hold the antenna upright (vertically) in the desired position, taking the following two things into account.
 1. If possible, position the antenna slightly higher than the device.
 2. Avoid placing it right next to metal objects, in a metal housing or near a frequency converter.
- Start up the device once to check the network strength. After start-up, the signal strength will appear on the display at the top. Wait a moment (1 to 2 minutes) before reading it, as the device needs time to search for the network.

3. Assessing signal strength

- If two or more bars are visible, the network connection is sufficient. The antenna can be mounted in this location. See the next page for instructions on how to do this.
- If fewer than two bars are visible: proceed to step 4.


4. Optimizing the antenna position

- Move the antenna in increments of no more than one metre. If no signal is detected at all, it is best to place the antenna on the outside of the stable.
- Adjust the angle of the antenna if necessary.
- Restart the device each time and check the network strength.

5. Repeat if necessary until two or more lines appear.

6.8 Additional points of interest

1. Drilling holes
 - Hold the antenna bracket in the previously determined position, then mark the two drill holes.
 - Use a drill with a 5 mm drill bit.
 - Drill the holes straight into the wall at the previously marked drill holes.
 - Remove drill dust from the holes (e.g. with a vacuum cleaner or by gently blowing into them).
2. Inserting plugs
 - Insert the 5 mm plugs into the drilled holes until they are flush with the wall.
3. Attaching the bracket with antenna
 - Hold the bracket in the correct position so that the holes in the bracket line up with the plugs.
 - Screw the screws through the bracket into the plugs.
 - Make sure the bracket is secure and straight.
 - Finally, apply tape to the connection between the antenna cable and the antenna (see image opposite). This increases reliability and extends the life of the connection.
4. Check
 - Check that the bracket is secure and safe.
 - Check the signal strength on the MS AutoHoofClean pump.

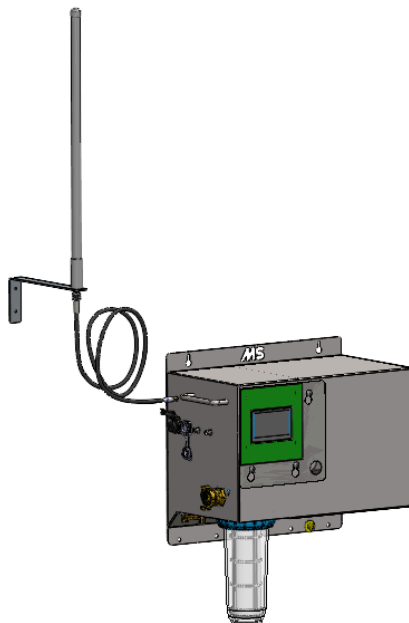6.9 Connecting power supply

The final step is to power up the system. To do this, connect the power cord to the AutoHoofClean pump. The power connection is located on the left side of the housing. Insert the connector into the socket and lock it by turning it so that the metal part is facing up. Check that the connector is securely locked by gently turning it.

Plug the system into a grounded outlet. Make sure the outlet is grounded and connected to a residual-current device (GFCI).


Once the plug is connected, the screen will turn on automatically.

7. Display and system settings

This chapter explains the AutoHoofClean pump display and how to set up the system. Please note: All images in this chapter are for illustrative purposes only and show the system with the language set to English.


Display overview and access to settings

When you are on the home screen, press the gear icon at the top right of the screen to open the settings menu.

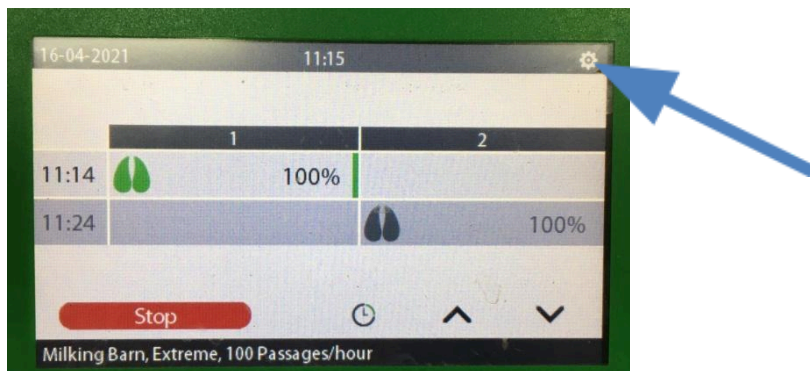This menu contains four options, which are briefly explained in the overview.

In this menu you can set, among other things:

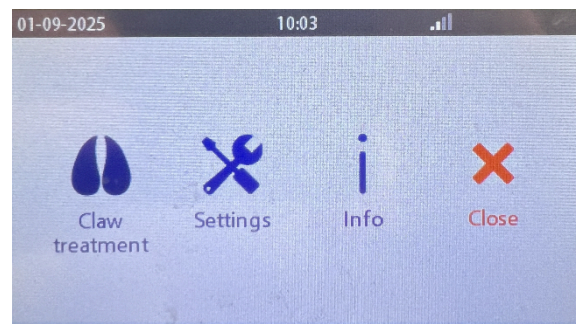
Claw treatment:

- The number of mats
- The start and end times of milking at a milking parlor
- The number of milkings per day at a milking parlor


By pressing the "More" button you can switch between the three screens below.

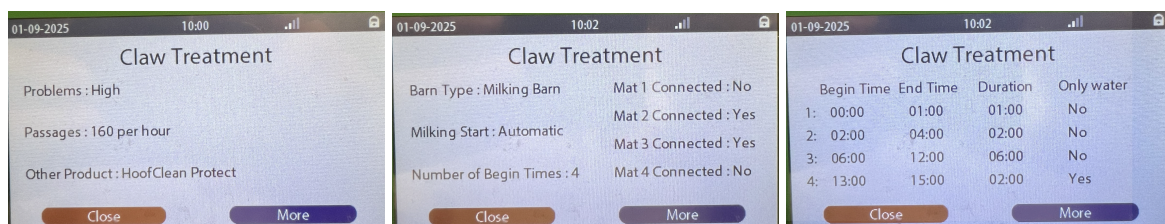
Settings:

This includes how to replace the filter, network, test, create a logbook and installer.

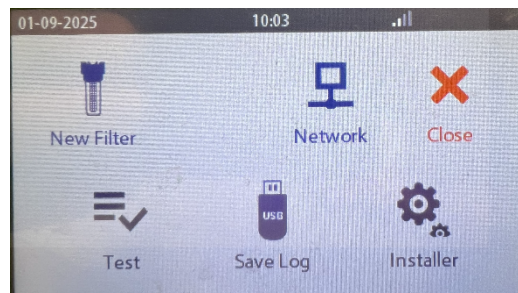Info:


Under info you will find general information about this product such as software and hardware versions.

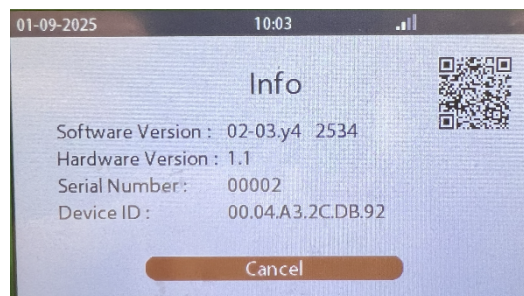Close:

Pressing "close" will take you back to the home screen.

To change settings or values, press the lock icon in the upper right corner for a few seconds to unlock it. Then enter your code [4444]. Once you're done, the adjustable values will turn blue, and you can make changes.


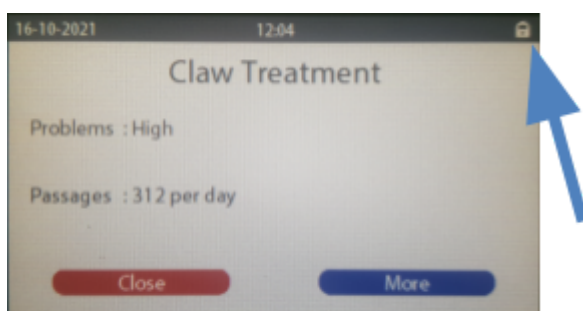
7.1 Setting language and time

It is important to first set the AutoHoofClean pump to the correct language and time before proceeding with other settings.

To do this, follow the steps below:

1. Press the gear icon at the top right of the screen.
2. Click on "Settings".
3. Select "Installer".
4. To unlock the settings, press and hold "Unlock" until a code field appears.
 - o Enter code 4444 to access the adjustments.
5. Adjust the values for "Language", "Day", "Month" and "Year" to your desired settings.
 - o You do this by clicking on the blue variables.
6. Change the values by clicking the left or right arrows.

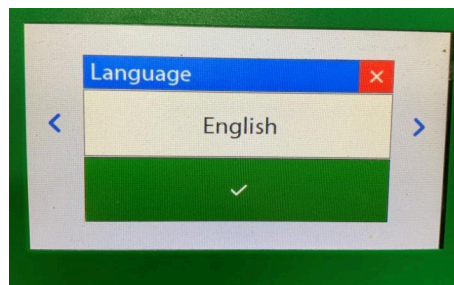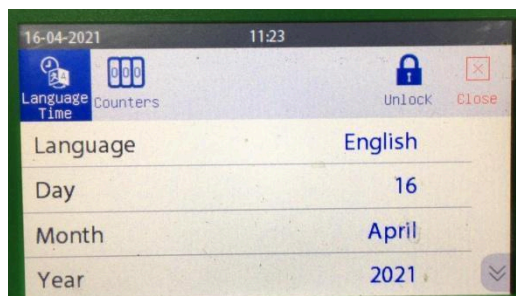8. Startup and testing


To use the system for the first time, the mats must be filled with water. This is done via the test screen.

The test screen can also be used to check whether fluid is actually coming out of a specific connection.

This chapter explains how to open and use the test screen.

-
1. Opening a test screen – step-by-step plan
 2. Press the gear icon at the top right of the screen.
 3. Go to the “Settings” menu.
 4. Press the “Test” icon at the bottom left.

The following screen will appear:

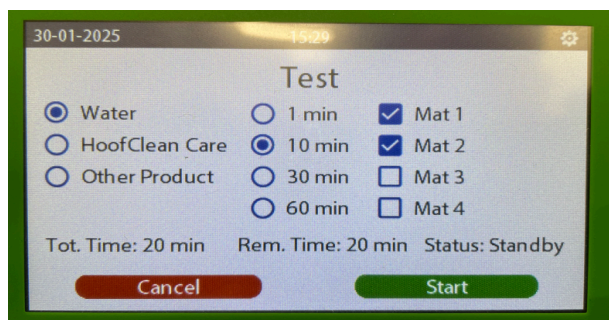You can test the system and fill the mats using the test screen above.

Setting up a test or filling program

1. Select the desired program

You can choose from the following three options:

- Water
- Water + MS HoofClean Care
- Water + Other product

2. Set the duration

Select the desired duration per mat:

- 1, 10, 30 or 60 minutes

Note: The time indicated is per mat. The system fills the mats in order, from mat 1 through mat 4.

3. Select the mats

Please indicate which mats you would like to test or fill.

Example:

If you choose the program *Water*, a period of time *30 minutes* If you select mats 1, 2, and 4, the system will fill these mats sequentially. The total time will then be 90 minutes (3 x 30 minutes).

Information on the screen

- Total Time: Shows the total duration of the test program.
 - Rem. Time: Shows the remaining time of the running program.
 - Status: Displays the current status of the testing process.
-

Recommendation for new mats

When a new mat is connected, it is recommended to fill it with the program *Water* for a duration of 30 minutes. Select the program for this *Water*, set the time to *30 minutes* and select the appropriate mats.

Do you want to check whether a product is being sucked in correctly?

Then choose a program with an added product (*Water + MS HoorClean Care* or *Water + Other product*) and observe the transparent suction hose. This allows you to see whether the correct product is being sucked in.

8.1 Select stable type

The AutoHoofClean pump is specially designed for use in two different types of barn: milking parlor or milking robot.

This chapter explains how to select and set the correct stable type.

Steps to select the stable type:

1. Press the gear icon at the top right of the screen to open the settings menu.
2. Select "Claw treatment".
3. Press and hold the lock icon in the top right corner of the screen for a few seconds.
4. Enter the code 4444 to access the adjustable parameters.

After entering the code, all adjustable values will turn blue.

5. Click "More" repeatedly until the stable type setting screen appears [see image].
6. Click on the blue value to select the desired stable type:
 - Milking parlor
 - Milking robot
7. Depending on the selected stable type, continue with the chapter:
 - Milking robot settings
 - Milking parlor settings

These chapters explain how to correctly configure the AutoHoofClean for your specific situation.

8.2 Milking robot settings

Note: If you are using the AutoHoofClean pump in a milking parlor, this chapter does not apply. In that case, proceed to the chapter: Milking Parlor Settings.

Navigate to the correct settings

Press the "More" button to scroll through the different settings until you reach the milking robot settings screen [see image].

Tip: Did you accidentally press "More" too many times? Keep clicking "More" until you reach the correct screen. The screens repeat in a loop.

Set number of mats

In the chapter "Connecting the mats to the AutoHoofClean pump," you connected the water hoses to the GeKa connectors on the bottom of the AutoHoofClean pump. Now you need to indicate which outlets are actually in use.


Activate or deactivate outputs

- For each output (mat) it must be set whether water and product should be supplied.
- Set the output to "Yes" if a mat is connected.
- Set the output to "No" if no mat is connected.

Method:

1. Press the blue highlighted text behind the relevant mat.
2. Select the desired setting: "Yes" or "No".
3. Confirm the choice by pressing the check mark.

Repeat this for all available exits (mats).

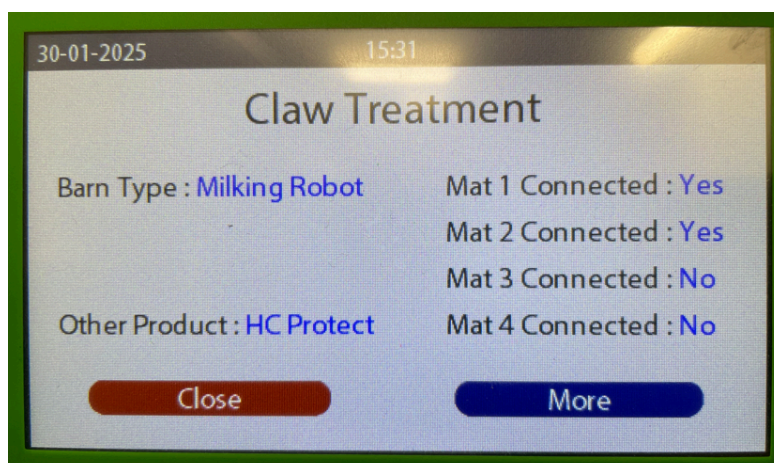8.3 Troubleshooting – milking robot

After setting the number of mats, the severity of the hoof problem must be set. This setting determines the intensity of the treatment via the AutoHoofClean pump.


Select problem level

1. Navigate to the correct screen by pressing “More” until you reach the “Problems” screen.
2. There are three adjustable levels:
 - Low
 - Normal
 - High
3. To adjust the value:
 - Press the blue highlighted value next to “Problems”.
 - Use the arrow keys to select the desired level.
 - Confirm your choice with the green check mark.
 - Or cancel the change with the red cross.

✓ Recommended starting setting:

Always set the problem to “High” when first using it.

As claw health improves over time, you can adjust the level to “Normal” or “Low” in consultation with your advisor.

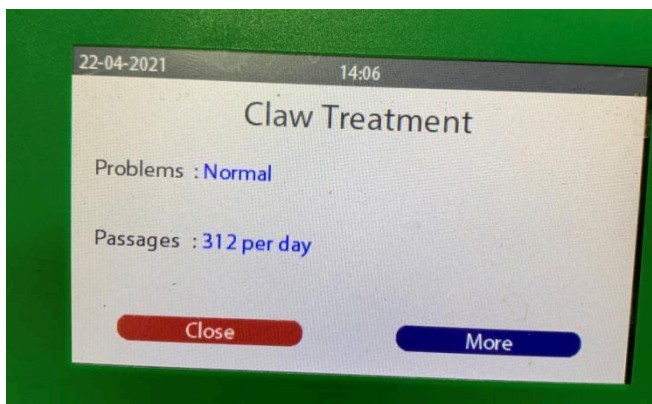8.4 Setting Passages

For the "Milking Robot" program to function correctly, it's important to set the correct number of passages per day. This information can be found in the milking robot's management package.

The setting concerns the number of cows that pass the milking robot per day.

Available settings:

1. Adjustable from 168 to 600 passages per day [with a step size of 48]
2. Always choose the value that is closest to the actual situation

Changing the number of passages:

1. Press the blue highlighted value to adjust it.
2. Use the left/right arrow keys to select the desired number of passages per day.
3. Press the green check mark to confirm the setting.

8.5 Set other product

To ensure the correct dosage of the disinfectant or protective agent used, it's important to select the correct product in the system. This prevents under- or overdosing and ensures the product works effectively.

Available product options:

You can choose from the following products:

3. Formades
4. MS HoofClean Protect
5. Product C

Setting up the product:

Setting the desired product is done via the following screen:


1. Navigate using "More" until you reach the screen with the "Other Product:" setting [see below].
2. Press the blue marked value behind "Other Product".
3. Select the correct product from the list.
4. Confirm your choice with the green check mark.

(You can cancel with the red cross.)

⚠ Ensure the selected product matches the product actually connected. Incorrect selection may result in incorrect dosing or damage to the mats/pump.

Click the blue value next to "Other Product." Then adjust it to the correct product using the left and right arrows. Finally, click the green check mark to save the setting. Once the product is selected, the AutoHoofClean pump knows which settings are required. You don't need to adjust anything else.

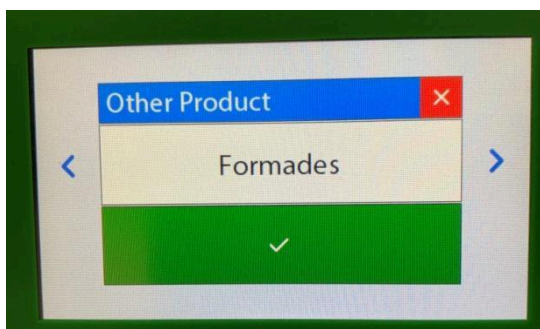9. Milking parlor settings

Please note: This chapter only applies if you are using the AutoHoofClean pump in a milking parlor.

If you use a milking robot, please refer to the chapter: Milking robot settings.

Navigate to the correct settings

1. Open the Claw Treatment menu.
2. Press "More" repeatedly until the correct settings screen appears.

Tip: Did you accidentally press "More" too many times? Press "More" again multiple times until you reach the correct screen – the screens will rotate in a loop.

Check or adjust stable type

Check that the correct barn type (Milking Parlour) is selected under "Barn Type".

To adjust this:

1. Press the lock icon for a few seconds to unlock the settings.
2. Enter the code 4444.
3. Click on the blue value behind "Barn Type" and select the correct type.
4. Confirm with the green check mark.

Set number of mats

In the chapter "Connecting the mats to the AutoHoofClean pump," you connected the water hoses to the GeKa connectors on the bottom of the AutoHoofClean pump. Now you'll set which outlets are actually in use.

Activate or deactivate outputs

- Set all outputs that have mats connected to them to "Yes".
- Set all unused outputs to "No".

Method:

1. Press the text "Yes" or "No" behind the relevant mat.
2. Change the setting to the desired value.
3. Confirm with the green check mark.

Repeat this step for each exit (mat) that is or is not in use.

NB! You can use a maximum of 2 mats in a milking parlor.

9.1 Starting milking – setting the function

When a milking parlor is selected as the barn type, the "Start Milking" function becomes available. This function determines when the AutoHoofClean pump should be active—i.e., when treatment takes place during milking.

You can choose between two modes:

- Automatic
- Manually

Setting the function

1. Go to the settings menu under Claw Treatment.
2. Find the "Start Milking" setting.
3. Press the blue value behind this setting.
4. Select Automatic or Manual depending on your preference.
5. Confirm your choice with the green check mark.

Continue to the appropriate paragraph based on your choice:

- Start automatic milking– for automatic activation of the AutoHoofClean pump during milking.
- Start manual milking– for manual operation of the AutoHoofClean pump during milking.

9.2 Start automatic milking

After selecting "Automatic" under Start Milking, press "More." A screen will then appear where you can set the following times, as shown here:

In this screen, you set the number of milking sessions per day. "Milk start times" determines how often you want to milk per day. You can choose from 1, 2, or 3 start times.

1. Duration: On the right you can set how long the AutoHoofClean pump should be active per milking moment.
2. Start times: Below you can set the start times of the milking sessions.
3. If you choose to set 3 milking times, a pop-up will appear asking whether only water should be used during the third time slot.
4. Click on the value in the pop-up and change it as desired.
5. Select "Yes" if you only want to use water during the third time slot.
6. Select "No" if you do not want this.

9.3 Starting manual milking

If you choose not to have milking take place at the same times every day, it is recommended that you start and stop the AutoHoofClean pump manually.

Setting up manual start

- Set the "Start Milking" value to "Manual".

Manual start and stop

- To start or stop the AutoHoofClean pump, perform a manual action via the control panel.
- Once milking is completed and you stop the AutoHoofClean pump, the system automatically rinses with water.

Note: It is possible to stop the rinse prematurely, but this is not recommended to ensure proper cleaning and functioning of the system.


9.4 Setting up the Claw Problem

After setting the mats, you need to determine the problem. This determines the severity of the claw problem and the intensity of the treatment.

1. Once you are in the screen below, press "More" to go to the issue settings screen.
2. You can choose from three adjustable levels:
 - Low
 - Normal
 - High
3. To adjust the setting:
 - Press the blue value next to "Problems".
 - Use the arrows to select the desired level.
 - Confirm your choice with the green check mark or cancel with the red cross.

Advice for first use: Always set the problem to "High".

If you notice that the severity of the problem is decreasing, consult with your advisor about lowering the level to "Normal" or "Low".

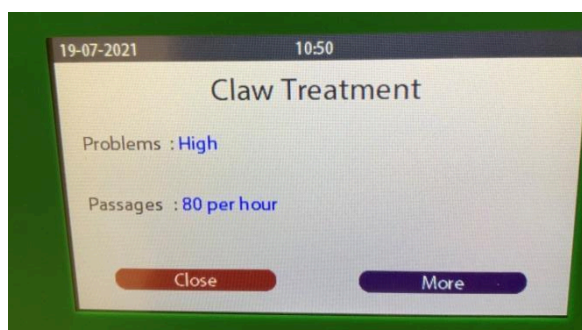9.5 Setting Passages

For the "Milking Parlor" program to function correctly, it's important to set the correct number of passages per hour. This information can be found in the farmer's records. The setting concerns the number of cows that pass through the milking parlor per hour. Available settings:

7. Adjustable from 40 to 200 passages per hour
8. In steps of 20 passages
9. Always choose the value that is closest to the actual situation

Changing the number of passages:

1. Press the blue highlighted value to adjust it.
2. Use the left/right arrow keys to select the desired number of passes per hour.
3. Press the green check mark to confirm the setting.

9.6 Set other product

For correct dosing of the disinfectant or protective agent, it is essential that the correct product is selected in the program.

You can choose from the following products:


- 10. Formades
- 11. MS HoofClean Protect
- 12. Product C

Select product

You can set the desired product in the following screen:

1. Navigate to the "Other Product:" screen [see below].
2. Press the blue value behind this option.
3. Select the product you want to use from the list.
4. Confirm your choice with the green check mark.


Ensure the product selected matches the product actually being used to ensure correct dosage and effectiveness.

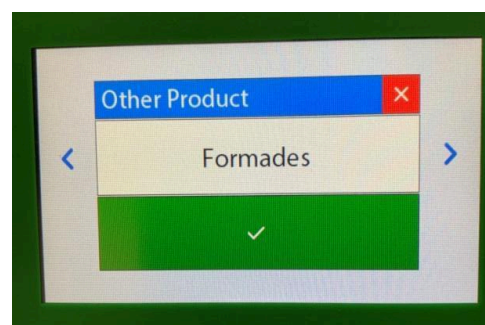
10. Meanings of icons on screen

It's important to know what the icons on the AutoHoofClean pump mean. The following images show some of the icons displayed when the AutoHoofClean pump is in use.


If two hooves are shown, this means that the system dispenses water + Hoofclean.

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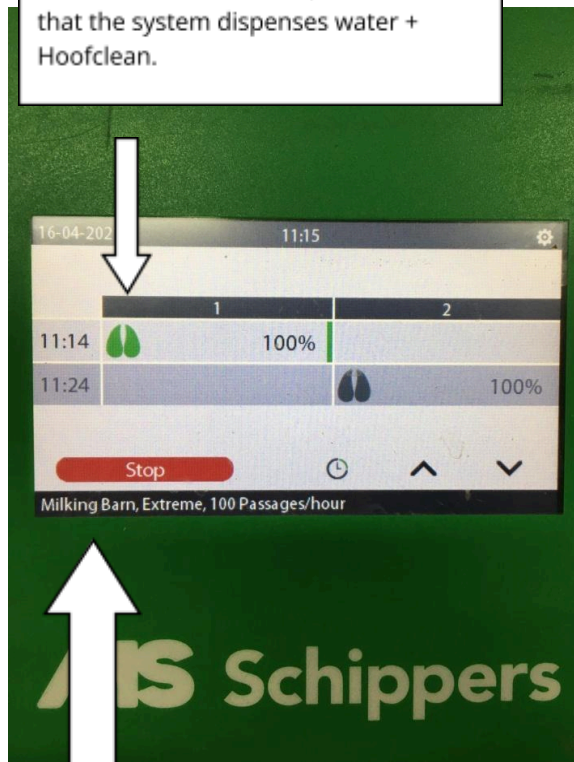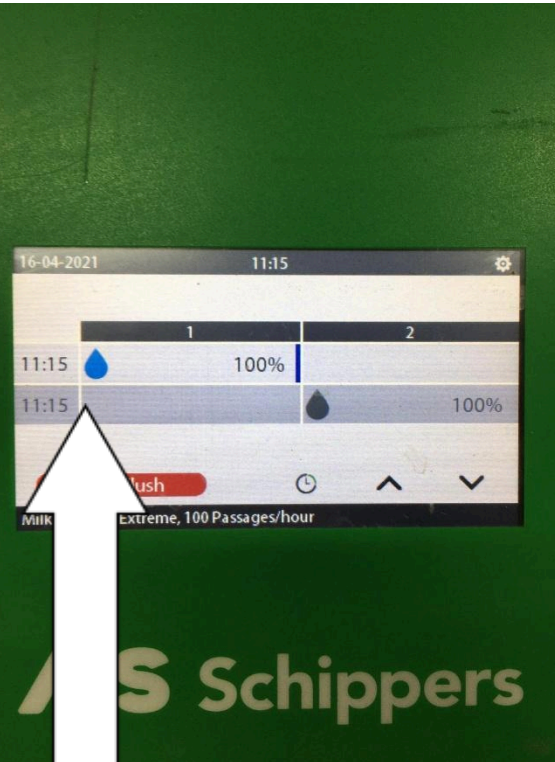If there is a water icon, this means that the system only provides water.

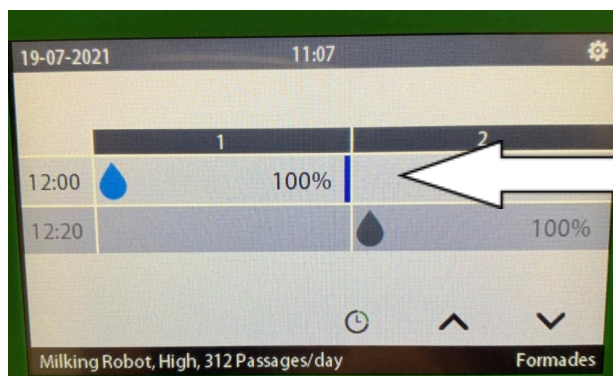If there is a bottle, this means that the system dispenses water + another product, for example MS HoofClean Protect.

This bar shows the progress of the cycle and gradually decreases as the cycle progresses.

At the bottom right, you can see which product you have selected. Ensure that this always corresponds to the product you are actually using.

11. Connect the AutoHoofClean pump via the HyCare app

1. Install app
Download the HyCare app from the App Store (iOS) or Google Play (Android) by searching for *HyCare*. Open the app, log in or create a new account, and follow the steps to select your location.
2. Turn on the pump
Plug the AutoHoofClean pump into the power outlet and wait until the light turns solid. Make sure your phone is nearby.
3. Connect pump
 - Open the Claw Health module in the app.
 - Choose AutoHoofClean pump.
 - Click the plus sign (+) to add a new device.
 - Follow the onboarding flow: scan the QR code on the pump or enter the pump ID manually.
 - Wait for the app to confirm that the pairing was successful.

The App is now ready to use.

12. MS HoofClean Care & other product empty notification

If your system detects that there is a (nearly) empty can, the following message will appear on the screen:

When you receive a message that the barrel is empty, you need to replace the product packaging. Follow the steps below carefully:

Steps to replace the barrel:

1. Unscrew the cap from the empty barrel.
Please note: leave the vacancy detector in the cap.
2. Place the new (full) barrel under the pump.
3. Remove the cap from the full barrel.
4. Screw the cap with the emptiness indicator onto the full barrel.
5. Confirm the replacement by pressing the "Barrel filled/replaced" button.

After this the system is ready for use again.

13. Read data via USB port

On the left side of the AutoHoofClean pump is a USB port. This port is used to read data from the system. The AutoHoofClean pump comes with a USB flash drive that can be used for this purpose.




You can read the following data via this USB port:

- Changes to the settings
- Consumption of products
- How the water meter works

Save log to USB

To save the log, follow the steps below:

1. Press the gear to go to the settings menu.
2. Go to "Settings".
3. Select "Save Log".
4. Insert the USB stick into the USB port on the left side of the AutoHoofClean pump.
5. Press "Save Log".
6. Wait until the system indicates that the USB stick can be safely removed.
7. Remove the USB stick.
8. You can read the saved data on a laptop or PC with programs such as Excel.

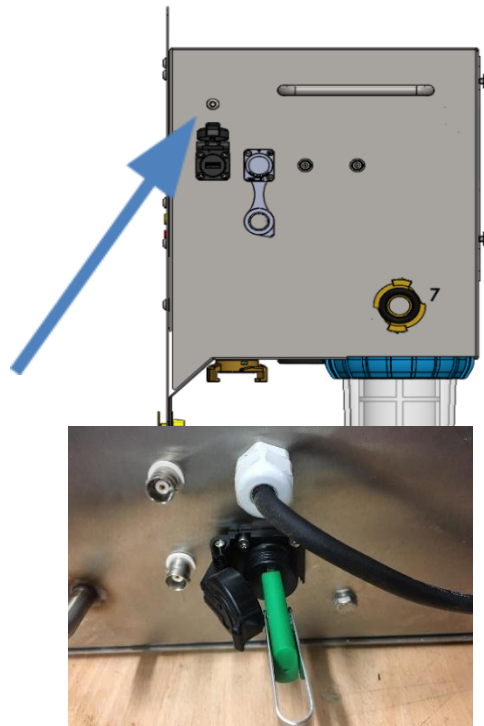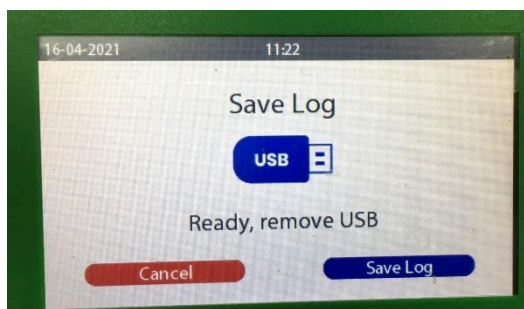
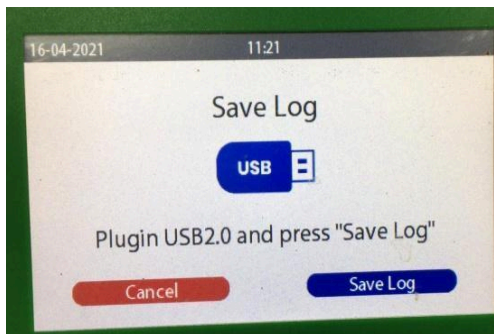
14. Screen installer

In the "Installer" screen you can set the date, time, language and screen brightness. There's also the option to read back the "Counters" in the program. Finally, under "Configuration," you can view and configure several settings, all of which are explained in this chapter.

14.1 Setting Language/Time

1. Go to "Settings" by pressing the gear at the top right.
2. Press "Installer" at the bottom right.
3. Press and hold the lock icon in the top right corner for a few seconds. Then enter the code 4444 to unlock additional features.
4. An icon for language and time will appear in the top left corner. You can make changes here.
5. Please note: All blue values can be changed without additional code.
6. Use the arrows on the right of the screen to scroll down and adjust further settings.

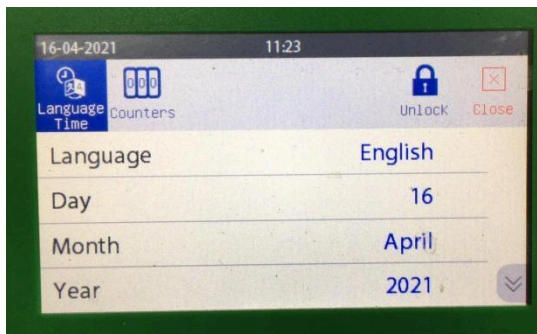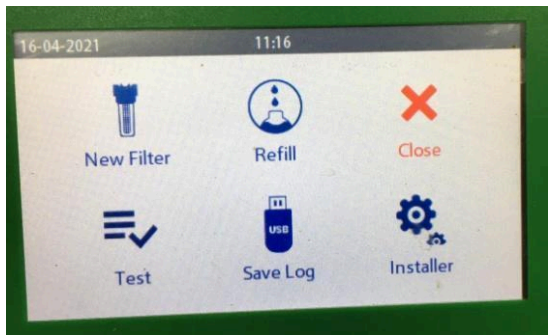14.2 Counter screen

The counter screen displays information about the number of pulses the AutoHoofClean pump has registered. This data is primarily informative and can help diagnose potential system malfunctions.



Here you see:

1. The number of pulses the system has received per day.
2. The total number of pulses since the system was started.
3. A counter that indicates whether the system has missed pulses. Missing pulses can indicate a possible leak.

Navigate to the counter screen


1. Tap "Counters" at the top of the screen.
2. Use the arrows to switch between the different screens.

Below you will find an overview of the relevant screens that are available.

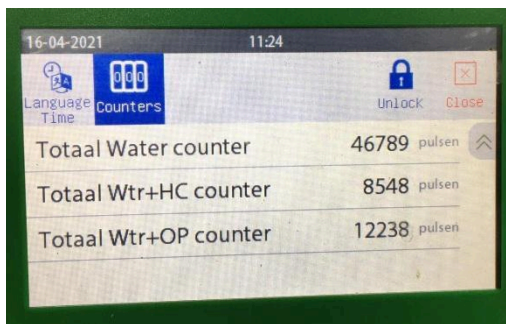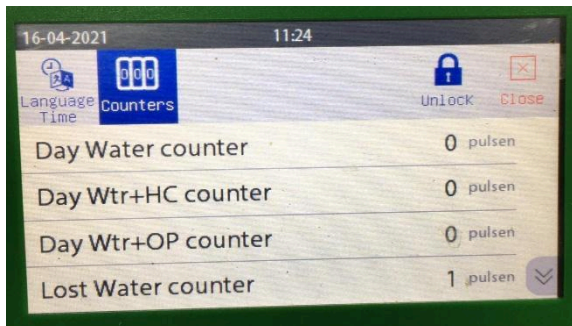14.3 Alarm level

1. Press and hold the lock icon at the top right and enter the code "4444" to access the advanced settings.
2. You can now view the alarm level and system configuration.
3. Use the arrows on the right side of the screen to navigate up and down through the menu.
4. In this menu you will find, among other things:
 - After how many liters the filter should be replaced.
 - How many liters are left before the filter needs to be replaced.

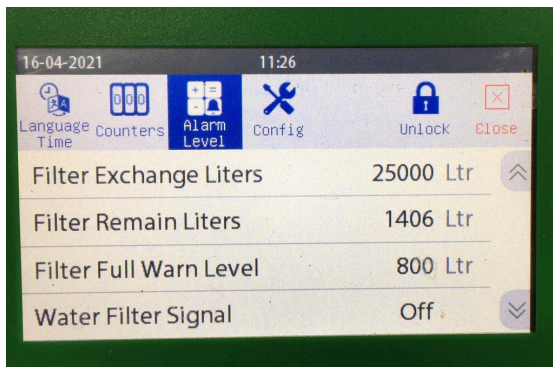15. Settings in the control panel (only applicable to milking parlors)

In the configuration screen two settings can be adjusted that only apply when using a milking parlor.

If the system is set up for use with a milking robot, you can skip this part of the configuration.

⚠ Please note: These settings are only relevant when manually starting and stopping the AutoHoofClean pump in a milking parlor.

Adjustable options:

1. Milk start code

This sets whether a code must be entered when starting the AutoHoofClean pump.

2. Milk stop code

This sets whether a code must be entered when the AutoHoofClean pump is stopped.


When you require a code, the user must enter the code "4444" when starting or stopping the system.

If you don't want to use a code, leave the values set to "No."

Access the control panel

1. Press the gear icon at the top right of the screen to open the settings menu
2. Select Installer

You will now enter the configuration screen (see image below).

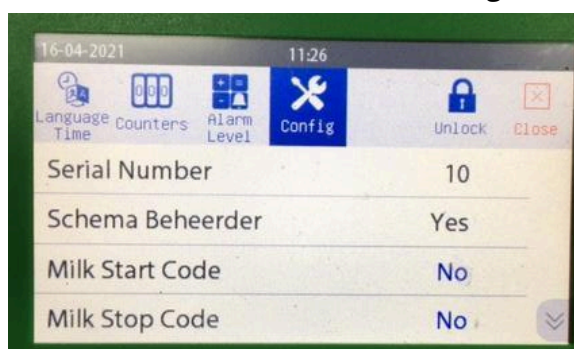3. Adjust the settings if necessary

Find the fields "Milk start code" and "Milk stop code".

- Choose whether to require a code when starting and/or stopping.
- Enter the desired setting (Yes/No).

16. View technical data via the "Info" menu



In the "Info" menu, you'll find important information about the AutoHoofClean pump. Here you can see, among other things:

- The software version of the pump
- The hardware version
- The serial number of the printed circuit board

This information is helpful to have on hand when you contact customer service. They may ask for it to help you more quickly.

Steps to find the information:

1. Press the gear icon at the top right of the screen to open the settings menu.
2. Select "Info".
3. The technical data screen appears.

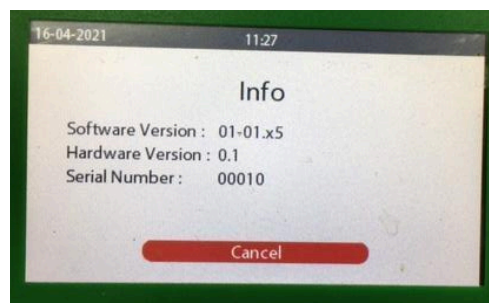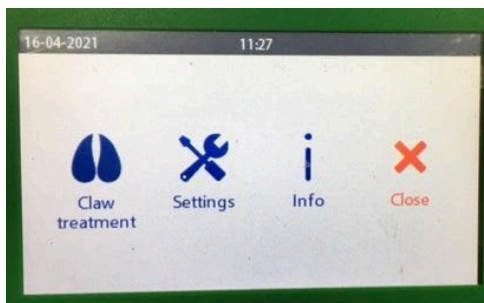
17. Maintenance of the Autohoofclean pump.

The AutoHoofCleanpmp is designed to be virtually maintenance-free. The only part that needs periodic replacement or cleaning is the water filter.

When the filter needs replacing or cleaning, the message "Filter needs cleaning" appears on the display. This message is displayed automatically after a preset number of liters of water has passed through the system.

The following section explains step by step how to properly replace or clean the filter.

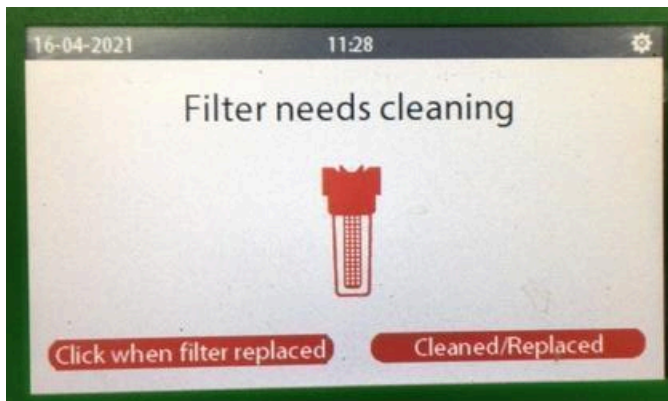18. Replacing the water filter (hardware)

The AutoHoofClean water filter is easily accessible at the bottom of the device. Follow the steps below to replace the filter safely and correctly.

⚠ Note: Always turn off the faucet before replacing the filter to prevent water damage. Without this step, water will continue to flow and the surrounding area may become wet.

Step-by-step plan:

1. Turn off the water supply
 - Turn off the tap to stop the flow of water to the appliance.
2. Take the pressure off the system
 - Run the test program briefly to reduce the pressure in the system.
3. Place a collecting bucket
 - Place a bucket under the filter to catch any residual water.
4. Unscrew the filter
 - Use the supplied filter wrench to unscrew the filter housing from the AutoHoofClean pump.
5. Remove the old filter cartridge
 - Remove the used filter cartridge from the filter holder.
6. Insert the new filter cartridge
 - Insert the new cartridge firmly into the filter holder.
7. Place the filter holder back on the appliance
 - Screw the filter holder back onto the AutoHoofClean-pump using the key.
8. Open the water supply
 - Turn on the tap to activate the flow of water to the appliance.
9. Check for leaks
 - Visually check that the filter is properly secured and that there are no water leaks.

18.1 Replacing the water filter (program)



After the water filter has been physically replaced, this action must also be registered in the AutoHoofClean pump software.

⚠ Important:

If the replacement is not registered, the system will continue to count and you may receive an early notification to replace the filter again – even though this is not necessary.

Step-by-step plan:

2. Open the settings menu
 - Press the gear icon at the top right of the screen.
3. Go to 'Settings'
 - In the menu, tap "Settings".
4. Open the filter menu
 - Select the "New filter" menu.
5. Confirm the replacement
 - Press the "Cleaned/Replaced" button to register that the filter has been replaced or cleaned.
6. Complete
 - The AutoHoofClean pump will now register the replacement and will not give an error message until the new filter actually needs replacing.

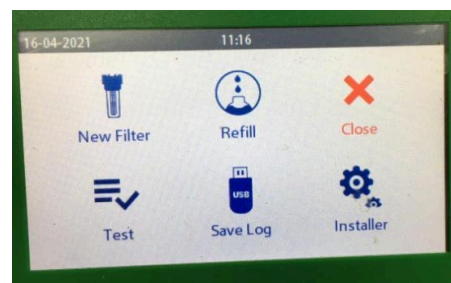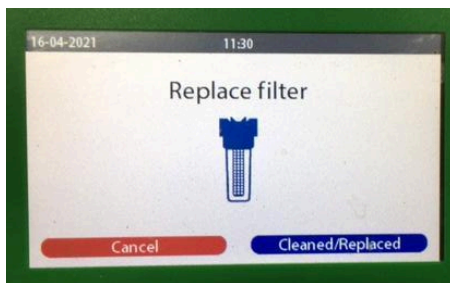
Completion of the installation

You have now successfully completed all the steps to properly use the AutoHoofClean pump.

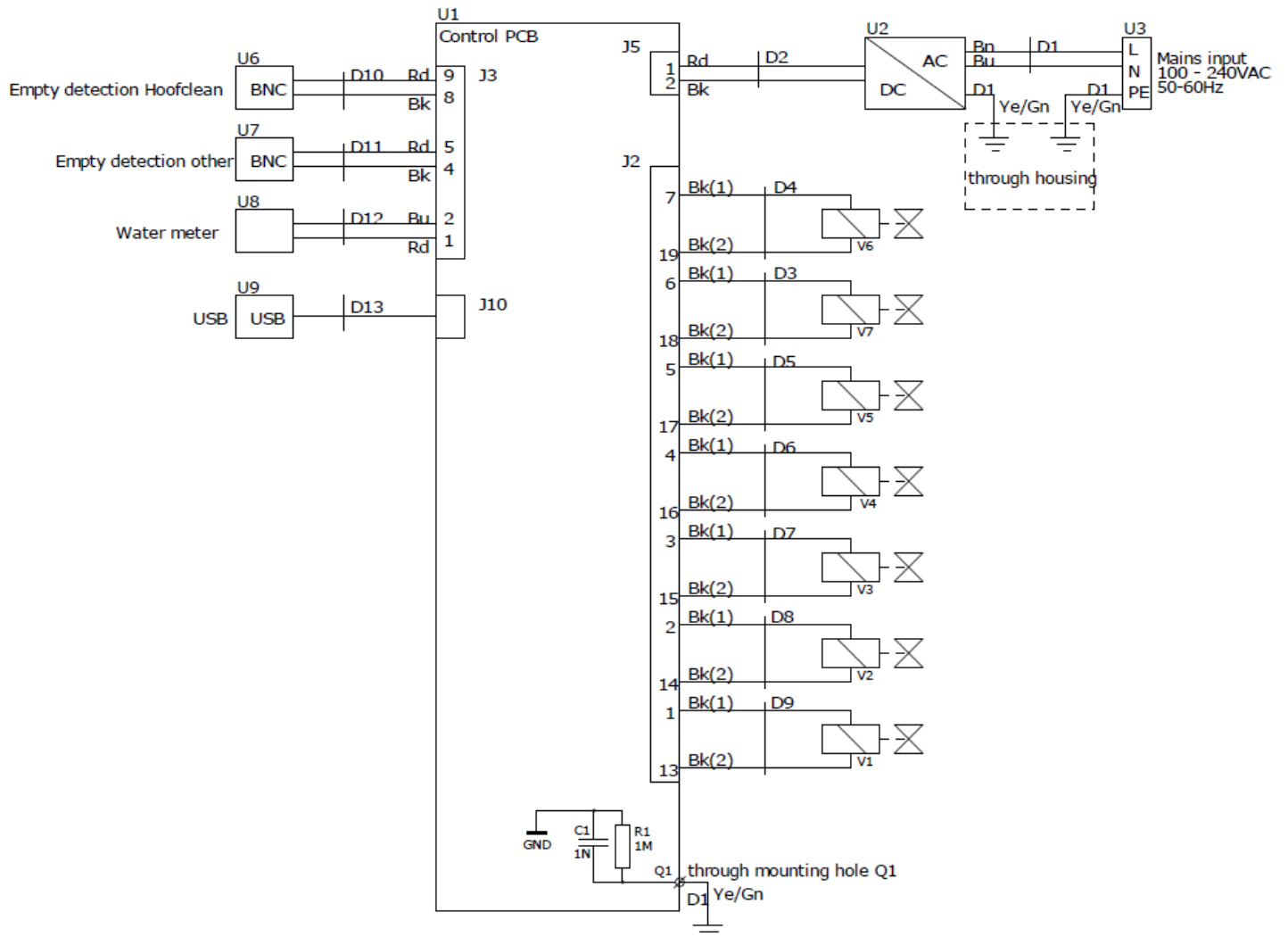We thank you for your confidence in our product. On behalf of The Schippers Group, we wish you ease of use and excellent results with your AutoHoofClean pump.

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© Copyright 2025